


# Triglycerides in Olive Oil - GC Application

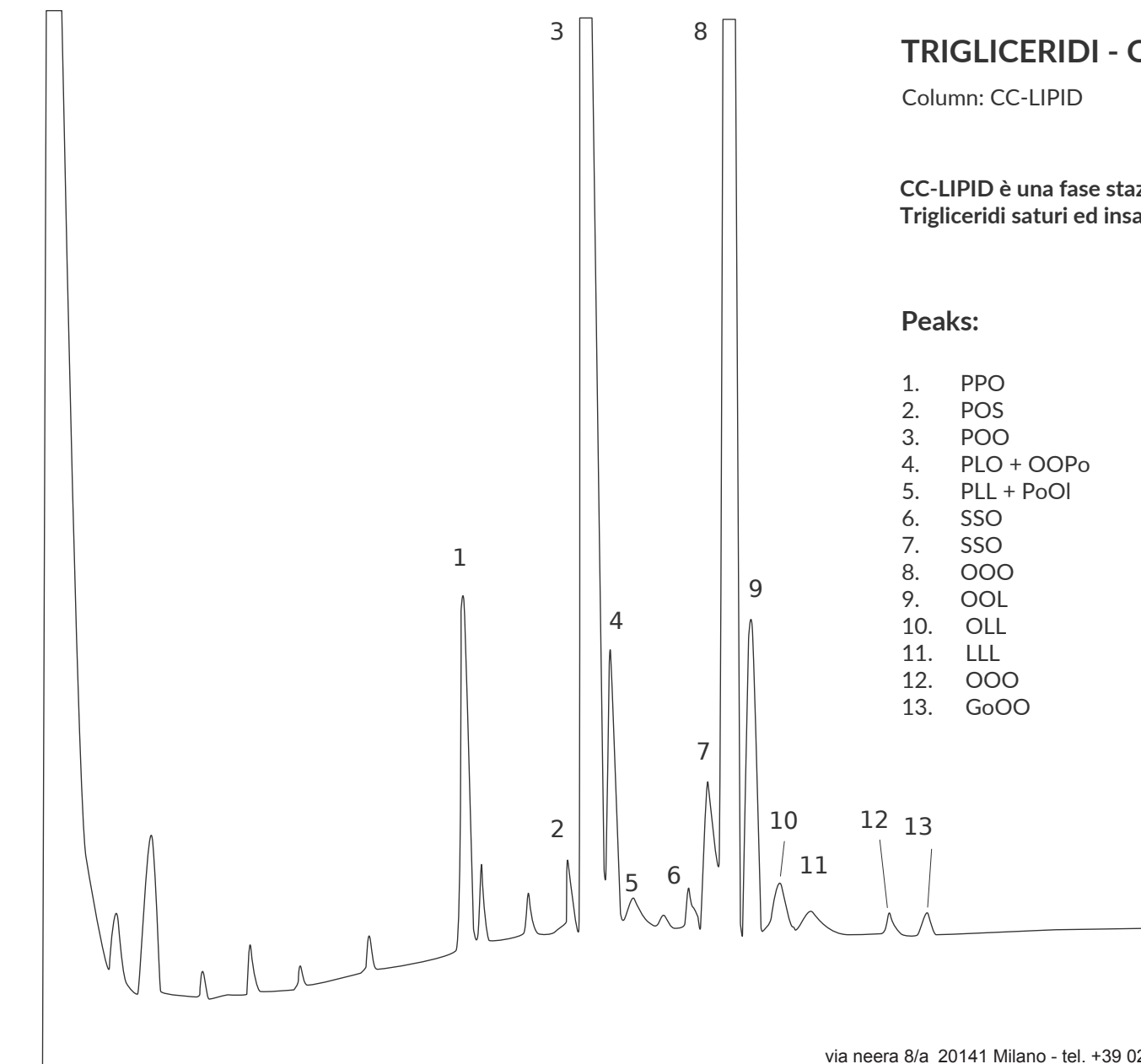

## TRIGLICERIDI - OLIO DI OLIVA

Column: CC-LIPID 25m x 0,32mm x 0,10µm

CC-LIPID è una fase stazionaria studiata per le analisi di Lipidi, Steroli e Trigliceridi saturi ed insaturi. Stabile fino a 370°C.

### Peaks:

1. PPO
2. POS
3. POO
4. PLO + OOPo
5. PLL + PoOI
6. SSO
7. SSO
8. OOO
9. OOL
10. OLL
11. LLL
12. OOO
13. GoOO

### Conditions:

Injection: Split 300°C, 1.0µL, 1:80 rapporto di Split

Detector: FID 300°C

Oven Temp.: 120°C, 5°C/min, 330°C

Carrier Gas: Hydrogen, 2mL/min